The Importance Of Pre-Cleaning Instruments At Point-Of-Use

Dr.Warda Abdelmoaty Azab
Treasurer of (ESSP)
HSPA Certified Professional (CRCST)

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Point-of-use cleaning is the first step in the decontamination process. It is performed at the point of use, where gross soil is removed from medical devices during and immediately after the procedure .


02

Point-of-Use Cleaning Or Pre-cleaning

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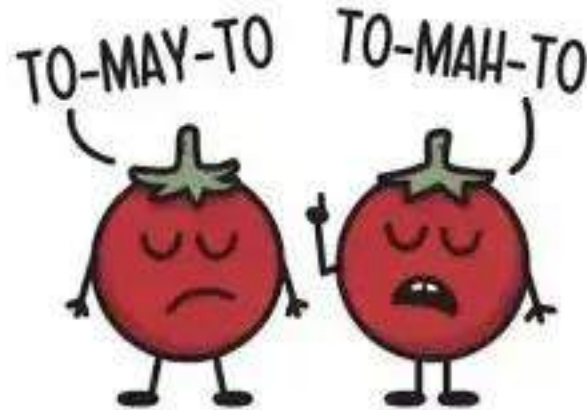- ✓ Blood
- ✓ Tissue
- ✓ Body fat
- ✓ Body Fluid

Why Remove Gross Soil

- ✓ Easier Cleaning
- ✓ Prevent corrosion
- ✓ Prevent dried debris
- ✓ Improves sterilization effectiveness

04 The importance of point-of-use cleaning:

- **Removal of gross soil** - Gross soil left on instruments not only makes them more difficult to clean but can damage instrumentation
- **Prevention of damage** - Soiled items must be treated in a manner that will prevent damage during transport to the decontamination area of the SPD.
- **Prevention of cross-contamination** - Soiled items must be safely transported from point-of-use to the decontamination area without cross-contaminating the environment between those areas.
- **Timeliness** - The prompt return of soiled or contaminated items to the SPD is essential for processing.
- **Keeping others safe**

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05 Reasons for point of use treatment


- Point-of-use treatment helps **prolong instrument life**.
- Dry soil and debris, especially in instruments with lumens and hard-to-reach crevices, are much more **difficult to remove than moist soil and debris**.
- When soil dries, instruments require more **aggressive cleaning methods** and may also take longer to clean. This can increase processing time and sometimes delay may impact instrument availability for use.
- Soil and excess moisture promote the formation of biofilm colonies. **Biofilm is highly resistant** to cleaning and disinfecting chemicals, so removing the causes of biofilm is essential.

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A group of microorganisms attached to surface and protected by a sticky layer

Why Is Biofilm Dangerous?

- ✓ Difficult to remove
- ✓ Resist to detergent
- ✓ Resist to disinfection
- ✓ Increase infection risk

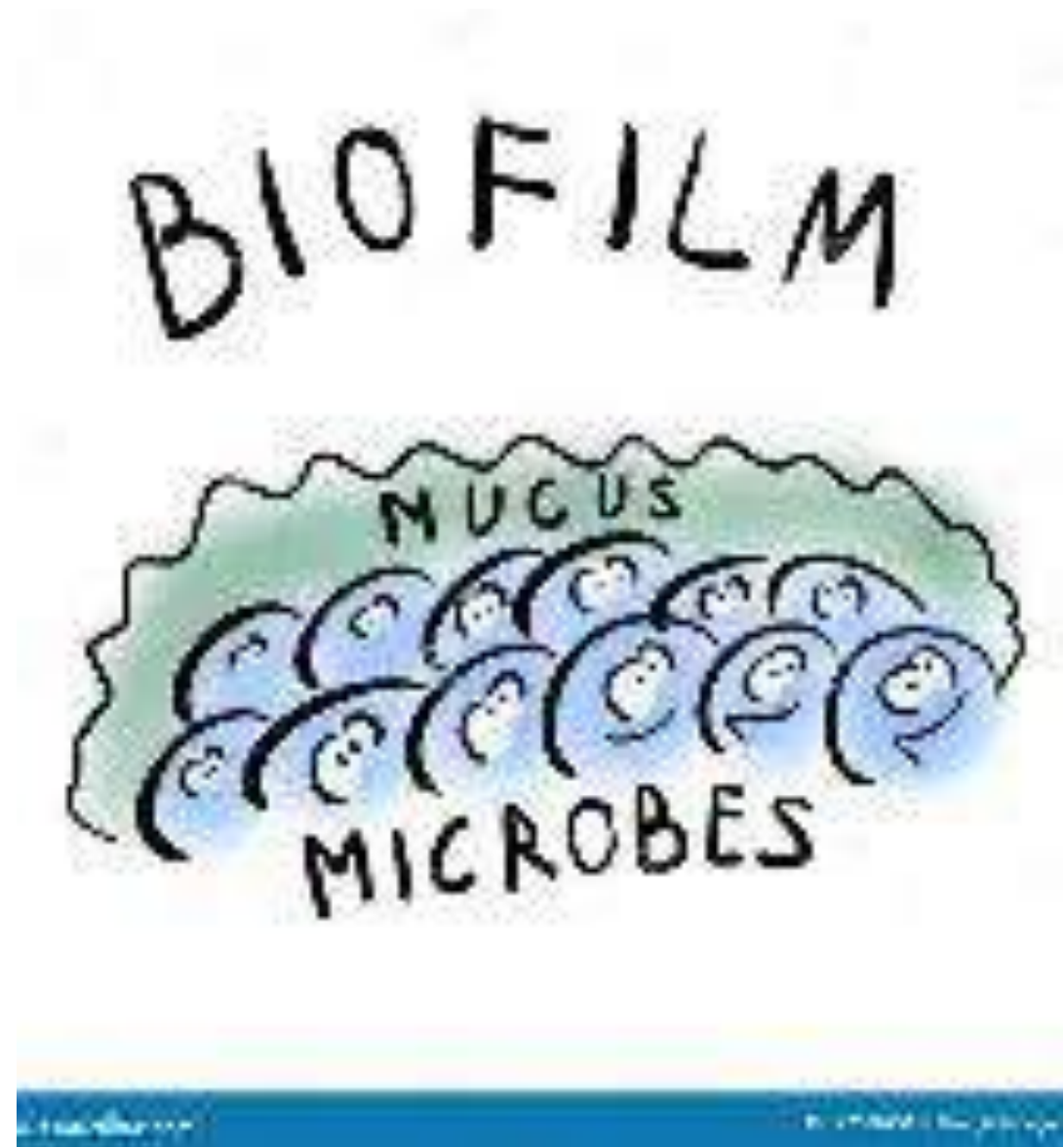Step-by-Step Point-Of-Use Cleaning

During and after procedure

- ✓ **Remove gross soil**
- ✓ **Keep instrument moist**
- ✓ **Open hinged instrument**
- ✓ **Separate sharp**
- ✓ **Remove disposables**
- ✓ **Keep sets together**

Handling sharps safely

Sharps includes

- ✓ **Scalpels**
- ✓ **Needles**
- ✓ **Scissors**
- ✓ **Bone instrument**

Rules

- ✓ **Place in puncture proof container**
- ✓ **Label as biohazard**
- ✓ **Avoid direct contact**
- ✓ **Never mix with other instrument**

Reusable Linen &Wast Separation

What to do?

- ✓ Remove reusable linen
- ✓ Place in correct bag
- ✓ Separate from instrument

Why?

- ✓ Prevent loss of small instrument
- ✓ Reduce contamination risk
- ✓ Improve sorting efficiency

Keeping Instrument Moist

Methods

- ✓ Moist towel
- ✓ Foam or gel spray
- ✓ Approved enzymatic product

DO NOT Over-Soak Instrument

Important warning

X Do not leave instrument in solution
for long time

Why?

- ✓ Can damage surface
- ✓ Causes corrosion
- ✓ Reduce instrument life

CJD Creutzfeldt-
Jakob Diseases
Precaution

If Suspected or Confirmed

- ✓ **Notify CSSD immediately**
- ✓ **Notify infection control**
- ✓ **follow hospital protocol for initial cleaning**

Repair & maintenance tagging

If instrument damage

- ✓ Tag for repair
- ✓ Separate from active sets
- ✓ Inform CSSD

Why?

- ✓ Prevent use the instrument
- ✓ Avoid procedures delays
- ✓ Maintain safety

Keeping Instrument Sets Together

Reules

- ✓ **Do not separate sets**
- ✓ **Keep all component together**
- ✓ **Arrange in order**

Why?

- ✓ **Prevent missing items**
- ✓ **Improve workflow**
- ✓ **Avoid delays in assembly**

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1. HSPA – Sterile Processing Technical Manual
2. CDC – Guideline for Disinfection and Sterilization in Healthcare Facilities
3. ANSI/AAMI ST79 – Comprehensive guide to steam sterilization

TOGETHER WE WIN